

3. Preparation of 0.05 (N) Mohr solution in 2 (N) H_2SO_4 (100 ml)Relative equivalent mass of $\text{FeSO}_4 \cdot (\text{NH}_4)_2\text{SO}_4 \cdot 6\text{H}_2\text{O} = 392.16$

Initial mass (g)	Final mass (g)	Mass of Mohr salt transferred (g)	Mass of Mohr salt to be taken (g)	Strength of Mohr solution
...	1.96	... (N)

[Dissolve the required amount of Mohr salt in 100 ml 2 (N) H_2SO_4 which is prepared by diluting 5.6 ml concentrated H_2SO_4 (36 N) to 100 ml carefully and then cooling.]

4. Titration and recording of e.m.f.: Volume of Mohr solution taken = 25 ml

No. of observations	No. of drops of $\text{K}_2\text{Cr}_2\text{O}_7$ (n) added	Observed e.m.f. (E_{cell}) volt or mV
1	0	
2	2	
3	4	
⋮	⋮	
⋮	⋮	
⋮	⋮	

5. Calibration of 1 cc of $\text{K}_2\text{Cr}_2\text{O}_7$ solution to number of drops : 1 cc = ... n' (drops say)

6. Graph plotting and Calculations

n	E_{cell} (volts)	ΔE_{cell} (volts)	Δn	$\Delta E_{\text{cell}}/\Delta n$ volt/drop
0 (n_1)	E_1	—	—	—
2 (n_2)	E_2	$E_2 - E_1$	$n_2 - n_1$	⋮
4 (n_3)	E_3	$E_3 - E_2$	$n_3 - n_2$	⋮
⋮	⋮	⋮	⋮	⋮
⋮	⋮	⋮	⋮	⋮

7. Calculation of strength of Mohr solution and formal potential of Fe (III)-Fe(II)—system

i. From the plot of E_{cell} vs. n , $E_{\text{Fe}^{3+}/\text{Fe}^{2+}}^0 = E_{\text{cal}} + (E_{\text{cell}})_{1/2} = \dots$ volt.

Expected value of $E_{\text{Fe}^{3+}/\text{Fe}^{2+}}^0 = +0.68$ volt in (M) H_2SO_4 medium


ii. From the plot of $\Delta E/\Delta n$ vs. n , number of drops of $\text{K}_2\text{Cr}_2\text{O}_7$ solution required = ... = ... ml

Hence strength of Mohr solution = ... (N).

Experiment 10

Potentiometric titration of Zn(II) solution by potassium ferrocyanide solution and also determination of the composition of Zn(II)-ferrocyanide complex.


Theory: The electrode potential for the redox system

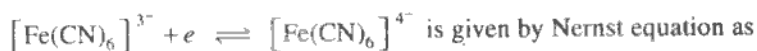

$$E = E^{\circ} + \frac{RT}{F} \ln \frac{C_{[\text{Fe}(\text{CN})_6]^{3-}}}{C_{[\text{Fe}(\text{CN})_6]^{4-}}}$$

where E° is the standard reduction potential of $[\text{Fe}(\text{CN})_6]^{3-} - [\text{Fe}(\text{CN})_6]^{4-}$ system.

If a clean platinum electrode is inserted into the solution of a mixture of $\text{K}_4[\text{Fe}(\text{CN})_6]$ and $\text{K}_3[\text{Fe}(\text{CN})_6]$, it will serve as an electrode. If we couple this electrode with a standard calomel (saturated) electrode, then the cell which is constructed for the experiment is

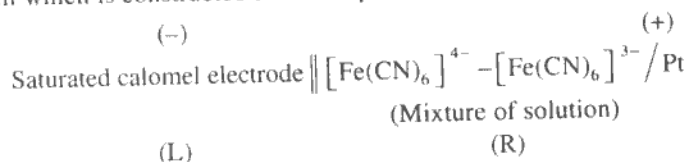

Pictorially it is represented as below:

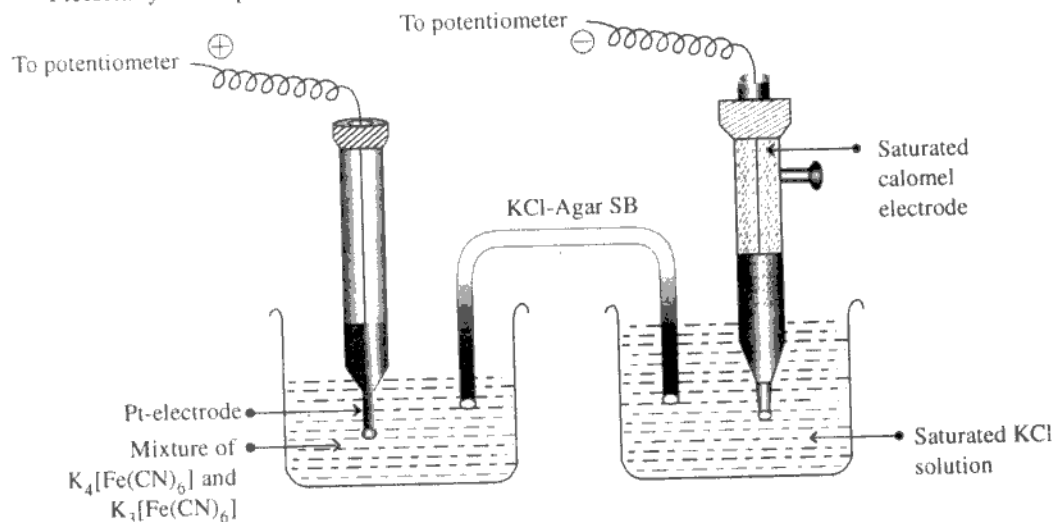


Fig. 6.25: Cell constructed for potentiometric titration of $\text{Zn}(\text{II})$ solution with $\text{K}_4\text{Fe}(\text{CN})_6$ solution

The observed e.m.f. of the cell is given by

$$\begin{aligned} E &= E_R - E_L \\ &= E^{\circ} + \frac{RT}{F} \ln \frac{C_{[\text{Fe}(\text{CN})_6]^{3-}}}{C_{[\text{Fe}(\text{CN})_6]^{4-}}} - E_{\text{cal}} = E^{\circ} - E_{\text{cal}} + \frac{RT}{F} \ln \frac{C_{[\text{Fe}(\text{CN})_6]^{3-}}}{C_{[\text{Fe}(\text{CN})_6]^{4-}}} \\ &= E^{\circ} - E_{\text{cal}} + 0.0591 \log \frac{C_{[\text{Fe}(\text{CN})_6]^{3-}}}{C_{[\text{Fe}(\text{CN})_6]^{4-}}} \quad [\text{at } 25^{\circ} \text{C}]. \end{aligned} \quad (1)$$

Here E_{cal} is the reduction potential of saturated calomel electrode.

When $\text{Zn}(\text{II})$ solution is added to the system, $[\text{Fe}(\text{CN})_6]^{4-}$ is removed and the proportions of $[\text{Fe}(\text{CN})_6]^{3-}$ increases. From equation (1) it is evident that with the progressive addition of $\text{Zn}(\text{II})$, the observed e.m.f. will gradually increase. At the equivalence point, there will be a sharp increase in e.m.f. due to sudden removal of all $[\text{Fe}(\text{CN})_6]^{4-}$.

If we plot $\Delta E/\Delta n$ vs. n (n is the number of drops of ZnSO_4 added) the nature of curve will be of the following type (Fig. 6.26).

From the equivalence point, corresponding to the peak, volume of ZnSO_4 required can be known. Volume of $\text{K}_4[\text{Fe}(\text{CN})_6]$ taken for estimation is also known. Therefore, strength of $\text{K}_4[\text{Fe}(\text{CN})_6]$ can be calculated as strength of ZnSO_4 is known (prepared).

Determination of composition of the complex formed

Let 10 ml $\text{K}_4[\text{Fe}(\text{CN})_6]$ is taken for titration.

Let the number of moles of $\text{K}_4[\text{Fe}(\text{CN})_6]$ present in 10 ml solution = a .

If at the equivalence point V ml ZnSO_4 solution is required, then V ml of ZnSO_4 solution combines with 10 ml $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution.

Let the number of moles of ZnSO_4 present in V ml solution = b .

Then find the ratio,
$$\frac{\text{No. of moles of } \text{ZnSO}_4}{\text{No. of moles of } \text{K}_4[\text{Fe}(\text{CN})_6]} = \frac{b}{a}$$

From this ratio we can conclude about the composition of the complex formed.

Materials

- 100 ml $\frac{M}{6}$ ZnSO_4 solution
- 100 ml $\frac{M}{100}$ $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution
- 100 ml $\frac{M}{100}$ $\text{K}_3[\text{Fe}(\text{CN})_6]$ solution
- Pt-electrode and saturated calomel electrode
- KCl-Agar salt bridge
- Potentiometer, (vii) Microburette, 100 ml beaker, pipette (10 ml), etc.

Procedure

- Prepare 100 ml $\frac{M}{6}$ ZnSO_4 solution, 100 ml $\frac{M}{100}$ $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution and 100 ml $\frac{M}{100}$ $\text{K}_3[\text{Fe}(\text{CN})_6]$ solution.
- Take 10 ml $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution and 10 ml $\text{K}_3[\text{Fe}(\text{CN})_6]$ solution (of order ~ 0.01 M) in a clean 100 ml beaker (use pipette) and insert the clean and dry Pt-electrode into the solution. Connect the saturated calomel electrode through KCl-Agar salt bridge as shown in Fig. 6.25. Then connect the cell to the terminals indicated in the potentiometer.
- Add one drop of ZnSO_4 solution from a microburette into the solution in the beaker—a white precipitate will be formed. Shake well and note the meter reading. Continue the titration adding 1 drop of ZnSO_4 solution at each time. Take 5 to 6 more readings after sharp change (equivalence point).

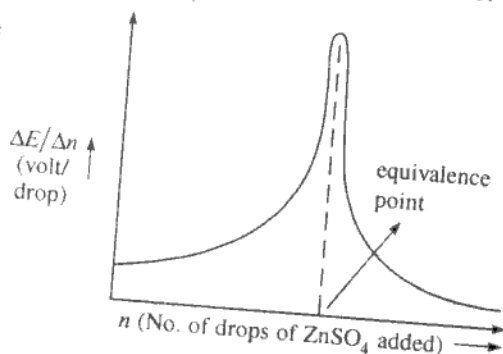


Fig. 6.26: Plot of $\Delta E/\Delta n$ vs. n (no. of drops of ZnSO_4 solution added for the titration of $\text{K}_4\text{Fe}(\text{CN})_6$ solution with ZnSO_4 solution).

- iv. Calculate $\Delta E/\Delta n$ and plot it against n . Detect the equivalence point in the graph and note the number of drops of ZnSO_4 required.
- v. Calibrate 20 drops of ZnSO_4 solution into volume (ml).
- vi. Calculate the strength of $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution. Also calculate the number of moles of ZnSO_4 and $\text{K}_4[\text{Fe}(\text{CN})_6]$. Then find ratio $\frac{b}{a}$ and conclude about the composition of $\text{Zn}(\text{II})$ ferrocyanide complex that is formed.

Results and Calculations

1. Recording of laboratory temperature

	Before experiment	After experiment	Mean
Temperature →	... °C	... °C	... °C

2. Preparation of 100 ml $\frac{M}{6}$ ZnSO_4 solution: Relative molar mass of $\text{ZnSO}_4 \cdot 7\text{H}_2\text{O} = 287.35$

Initial mass (g)	Final mass (g)	Mass of ZnSO_4 transferred (g)	Mass of ZnSO_4 to be taken (g)	Strength of ZnSO_4 solution
...	4.789	... (M)

3. Preparation of 100 ml 0.01 (M) $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution: Relative molar mass of $\text{K}_4[\text{Fe}(\text{CN})_6] \cdot 3\text{H}_2\text{O} = 422.26$.

Initial mass (g)	Final mass (g)	Mass of $\text{K}_4[\text{Fe}(\text{CN})_6]$ transferred (g)	Mass of $\text{K}_4[\text{Fe}(\text{CN})_6]$ to be taken (g)	Strength of $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution
...	0.4222	... (M)

4. Preparation of 100 ml 0.01 (M) $\text{K}_3[\text{Fe}(\text{CN})_6]$ solution: Relative molar mass of $\text{K}_3[\text{Fe}(\text{CN})_6] = 329.16$

Initial mass (g)	Final mass (g)	Mass of $\text{K}_3[\text{Fe}(\text{CN})_6]$ transferred (g)	Mass of $\text{K}_3[\text{Fe}(\text{CN})_6]$ to be taken (g)	Strength of $\text{K}_3[\text{Fe}(\text{CN})_6]$ solution
...	0.3291	... (M)

5. Titration and recording of e.m.f.'s

Volume of $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution taken = 10 ml
 Volume of $\text{K}_3[\text{Fe}(\text{CN})_6]$ solution taken = 10 ml } in a 100 ml beaker

Observation number	No. of drops of ZnSO_4 solution added	Observed e.m.f. (in V or mV)
1	1	...
2	2	...
3	3	...
4	4	...
⋮	⋮	...

6. Calibration of number of drops of ZnSO_4 solution to volume (ml)

Observation No.	No. of drops of ZnSO_4 solution (from microburette)	Volume of 20 drops of ZnSO_4 solution (ml)	Mean volume of 20 drops (ml) ZnSO_4 solution	No. of drops for 1 ml solution
1	20	...	V ml (say)	$\frac{20}{V} = n'$ (say)
2	20	...		
3	20	...		

7. Graph plotting and calculation of strength of $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution

No. of drops of ZnSO_4 solution added	Observed e.m.f. E (volt)	ΔE (volt)	Δn	$\frac{\Delta E}{\Delta n}$ volt/drop
1 (n_1)	E_1	—	—	—
2 (n_2)	E_2	$E_2 - E_1$	$n_2 - n_1$...
3 (n_3)	E_3	$E_3 - E_2$	$n_3 - n_2$...
...
...

From the graph, obtained by plotting $\Delta E/\Delta n$ vs. n , it is found that the number of drops of ZnSO_4 solution, at equivalence point (corresponding to the peak) = n'' (say). So, the volume of ZnSO_4 required for 10 ml $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution = $\frac{n''}{n'} \text{ ml} = x \text{ ml}$ (say).

Therefore, the strength of $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution

$$= \frac{x \times S}{10} (\text{M}) = y (\text{M}) \text{ (say), where } S \text{ is the strength of } \text{ZnSO}_4 \text{ solution in molarity.}$$

Compare this strength, obtained potentiometrically, with the strength of $\text{K}_4[\text{Fe}(\text{CN})_6]$ solution you have prepared by accurate weighing.

Determination of Composition of $\text{Zn}(\text{II}) - \text{Fe}(\text{CN})_6^{4-}$ complex

$$\text{Number of moles of } \text{ZnSO}_4 = \frac{x \times \text{molar concentration of } \text{ZnSO}_4 \text{ solution}}{1000} = b$$


Again, number of moles of $\text{K}_4[\text{Fe}(\text{CN})_6]$ that combine with ZnSO_4

$$= \frac{10 \times \text{molar concentration of } \text{K}_4[\text{Fe}(\text{CN})_6]}{1000} = \frac{y}{100} = a.$$

$$\text{Thus, the ratio} = \frac{\text{No. of moles of } \text{ZnSO}_4}{\text{No. of moles of } \text{K}_4[\text{Fe}(\text{CN})_6]} = \frac{b}{a}.$$

Expected value of the ratio is $\frac{3}{2}$.

Hence, the precipitation reaction that takes place in neutral medium according to the following equation,

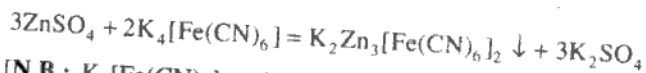

[N.B.: $\text{K}_4[\text{Fe}(\text{CN})_6]$ and $\text{K}_3[\text{Fe}(\text{CN})_6]$ must be freshly prepared. Otherwise their concentrations will change in presence of air and the result will be inaccurate].